



Put yourself in leading clinician's shoes

High immediate bond strength allows me the confidence to perform final adjustments without disrupting the bonded interface and unless you want to wait 24 hours to adjust, finish or polish you need - Multilink.

Dan Poticny, DDS
Dallas, TX

I saw it for myself - an additional 15 seconds of primer application provides five times the bond strength of competitive materials - it's worth it to secure my patient's restorations.

James Simon DDS
Memphis, Tennessee

Proven technologies and techniques are a must! With Multilink Automix, not only do I have proven technology, I have it in a faster more convenient technique allowing me the opportunity to bond restorations as beautifully as ever, but in a far more convenient fashion.

Robert Conrad, DDS
Houston, TX

I use Multilink Automix with my CAD/CAM restorations. Multilink's universal indications for both all-ceramic and zirconia restorations simplifies my workflow and the high bond strength gives me confidence that my restorations will be clinically successful for years to come.

Alex Touchstone DDS
Hattiesburg, MS

**Multilink[®]
Automix**
UNIVERSAL ADHESIVE CEMENT


**ivoclar[®]
vivadent[®]**

Strong


Multilink Automix provides high immediate bond strengths as well as long term security to meet you and your patient's expectations. Multilink Automix contains hydrolytically stable phosphonic acids that provide not only high bond strengths, but results that will last. Strength matters, choose Multilink Automix.

Shear Bond Strength Test Results


All-ceramic

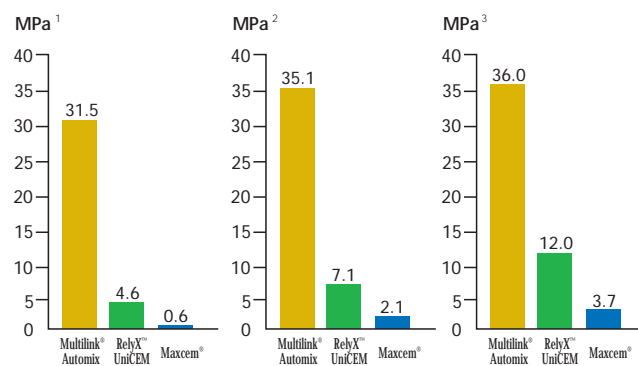

Zirconia

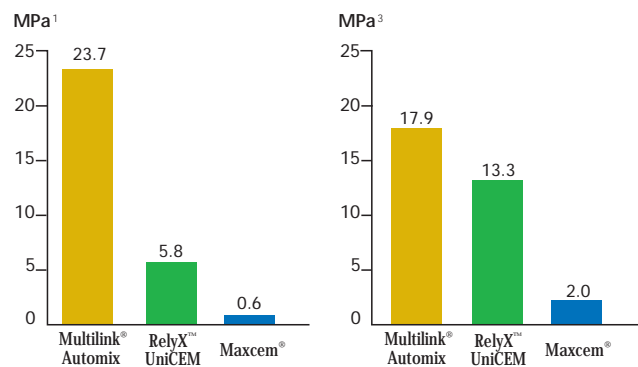

Metal-based[†]

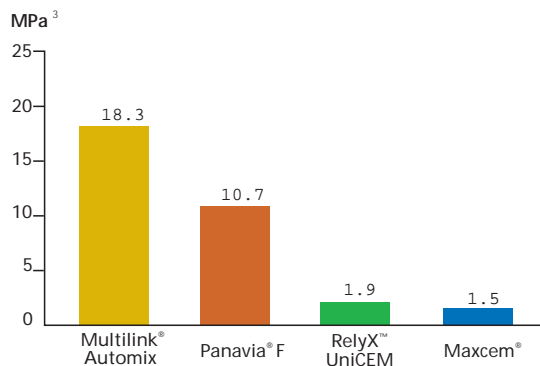


Fast

Multilink Automix cement is a universal adhesive cement that simplifies the cementation process for you.

Multilink Primer is the self-etching, self-curing primer used in conjunction with Multilink Automix cement. A quick application (15 seconds) seals the dentin, provides good marginal adaptation, and high bond strength.



Multilink Automix is available in three shades to meet the clinical indication and esthetic demands.

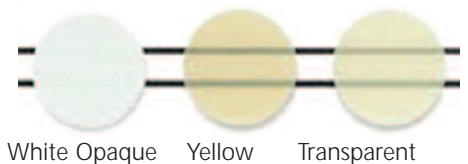

Applied per manufacturers instructions; Light-cured with high output quartz halogen curing light (>1,000 mW/cm²); Full report available upon request

¹Ivoclar Vivadent Applied Testing Center
Amherst, NY, February 2005

²James Simon, DDS, Med
Waldemar deRijk, DDS, Ph.D.
Director of Esthetic Dentistry and Clinical Research
University of Tennessee, February 2005

³Dr. Carlos Munoz
Chairman of Restorative Dentistry
University at Buffalo, February 2005

[†]10 minute bond strength results; self-cure; partial light-cure for 10 seconds per surface

Multilink Automix is a trademark of Ivoclar Vivadent, Panavia F is a trademark of Kuraray Co. Ltd, RelyX Unicem is a trademark of 3M ESPE, Maxcem is a trademark of Kerr Corporation.

Easy

Introducing... The *Ultimate*
Universal Resin Cement



Before



prime



seat



clean up



cure



After

Universal

Multilink Automix is a truly multi-purpose adhesive cement, indicated for cementation of indirect restorations (inlays, onlays, crowns, bridges, and posts) made of:

- All-ceramic
- Zirconia and aluminum oxide
- Metal and metal ceramics (PFM)
- Fiber reinforced composite



Delivery Forms



Multilink Automix System Packs
591406 Yellow
591405 Transparent
591407 White Opaque

Contents:

(1) – Automix Syringe (9g)
(1) – Primer A (3g)
(1) – Primer B (3g)
(1) – Metal/Zirconia Primer (5g)
(1) – Monobond-S Silane Coupling Agent (5ml)
(15) – Automix Tips
Various Accessories (brushes, mixing pad, applicator tips)



Multilink Automix Refills
591409 Yellow
591408 Transparent
591410 White Opaque

Contents:

(1) – Automix Syringe (9g)
(15) – Automix Tips



576825 Multilink Primer A/B Refill
(1) – Primer A (3g)
(1) – Primer B (3g)



579225 Multilink Primer A/B Mixing Pad



532888 Monobond-S™ Refill (5ml)

576827 Multilink Metal/Zirconia Primer

592225 Multilink Automix Tip Refill
(15) – Automix Tips

592225 Multilink Applicator Refill
(50) – Brushes

603018 Multilink Root Canal Tip Refill (Available January 2007)
(5) – Applicator Tips



USA

175 Pineview Drive
Amherst, NY 14228
716.691.0010
1.800.533.6825
716.691.2285 Fax
www.ivoclarvivadent.us
mail@ivoclarvivadent.us

CANADA

2785 Skymark Avenue, Unit 1
Mississauga, Ontario L4W 4Y3
905.238.5700
1.800.263.8182
905.238.5711 Fax
www.ivoclarvivadent.ca
mail@ivoclarvivadent.ca

100% CUSTOMER SATISFACTION
GUARANTEED!

©2006 Ivoclar Vivadent. Call us toll free at 1-800-533-6825 in the U.S., 1-800-263-8182 in Canada.
Multilink is a registered trademark of Ivoclar Vivadent.
Clinical photography courtesy of Dr. Robert Ritter, Palm Beach Gardens, Florida

594992 10/2006

Multilink®
Automix
UNIVERSAL ADHESIVE CEMENT

ivoclar
vivadent®